FESTO

Robotino®
Handbuch

Robotino®
Manual

mit CD ROM
CD-ROM included

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This mobile robot system Robotino® has been developed and produced solely for vocational and further training purposes in the field of automation and technology. The company undertaking the training and/or the instructors is/are to ensure that trainees observe the safety precautions described in the manuals provided. Festo Didactic herewith excludes any liability for damage or injury caused to trainees, the training company and/or any third party, which may occur if the system is in use for purposes other than purely for training, unless the said damage/injury has been caused by Festo Didactic deliberately or through gross negligence.

Bestell-Nr./Order No.: 544305
Stand/Status: 03/2007
Autor/Author: Ralph-Christoph Weber
Grafik/Graphics: Doris Schwarzenberger

© Festo Didactic GmbH & Co. KG, 73770 Denkendorf, Germany, 2007
Internet: www.festo-didactic.com
e-mail: did@de.festo.com

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The copying, distribution and utilisation of this document as well as the communication of its contents to others without express authorisation is prohibited. Offenders will be held liable for the payment of damages. All rights reserved, in particular the right to carry out patent, utility model or ornamental design registration.

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1.	Introduction	49
1.1	Training contents	49
1.2	Learning aims	49
1.3	Important notes	50
1.4	Duty of the operating authority	51
1.5	Duty of trainees	51
1.6	Risks involved in dealing with Robotino®	51
1.7	Warranty and liability	52
1.8	Intended use	52
2.	Notes on safety	53
3.	Technical data	54
4.	Transport/Unpacking/Scope of delivery	55
5.	Design and function	56
5.1	Robotino®	56
5.2	Chassis and command bridge	57
5.3	The drive unit module	60
5.4	The camera module	62
5.5	The controller unit	63
5.6	I/O circuit board module	63
5.7	Power supply/battery charger	64
5.8	Sensors	64
5.9	Membrane keypad and display	68
5.10	Wireless LAN access point	69
5.11	The compact flash card	69
5.12	The I/O interface	70
6.	Commissioning	71
6.1	Initial steps	71
6.2	Switching On and Off	71
6.3	The display functions	72
6.4	Testing the demo programs	73
6.5	Executing demo programs on solid ground	75
6.6	Controlling the Robotino® with Robotino®View	75
6.7	Setting up a WLAN connection	75
6.8	Testing the WLAN connection	76
6.9	Working with several Robotinos®	81
6.10	Charging the batteries	83
6.11	Replacing the batteries	85
6.12	Installing additional sensors	88
6.13	Connecting the sensors	89
6.14	C++ libraries	90
7.	Content of the CD-ROM	92
7.1	Operating instructions and Data sheets	92

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1. Introduction

The Festo Didactic Learning System for Automation and Technology is designed to meet a number of different training and vocational requirements. The mobile robot system Robotino® facilitates industry-orientated vocational and further training and the hardware consists of didactically suitable industrial components.

The mobile robot system Robotino® provides you with an appropriate system for practice-orientated tuition of the following key qualifications

- Social competence,
- Technical competence and
- Methodological competence

1.1

Training contents

Training contents covering the following subjects can be taught:

- Mechanics
 - Mechanical construction of a mobile robot system
- Electrical
 - Control of drive units
 - Correct wiring of electrical components
- Sensors
 - Sensor-guided path control
 - Collision-free path control with distance sensors
 - Path control by image processing of webcam pictures
- Feedback control systems
 - Control of omnidirectional drives
- Use of communication interfaces
 - Wireless LAN
- Commissioning
 - Commissioning of a mobile robot system

1.2

Learning aims

Beside industrial robot technology grows the market and therefore the importance of mobile robot systems and service robots. With Robotino® education responds to this technical and economic trend.

With Robotino® the learning aims below can be achieved.

The trainees

- learn to manage an electrically controlled drive unit
- know the fundamentals, the construction, the measurement of values and the parameterisation of a direct current motor control
- know the fundamentals of electrical drive technology
- understand an omnidirectional 3-axis drive and know to operate it
- are able to realise the commissioning of a mobile robot system using the Robotino® as an example
- can move the Robotino® in different directions
- can realise a sensor-guided path control
- can integrate image processing into the control of the Robotino®
- can develop a sensor-guided autonomous path control of the Robotino® using object recognition for example by colour.

Furthermore the following additional training aims can be achieved

- mounting and integration of additional sensors into the mobile robot system Robotino®
- can mount and integrate handling devices into the system
- are able to program their own navigation and control algorithms using a C++ library
- create a autonomous navigation of the Robotino®.

1.3

Important notes

The basic requirement for safe use and trouble-free operation of the mobile robot system Robotino® is to observe the fundamental safety recommendations and regulations.

This manual contains important notes concerning the safe operation of the mobile robot system Robotino®.

The safety recommendations in particular must be observed by anyone working on the mobile robot system Robotino®.

Furthermore, the rules and regulations for the prevention of accidents applicable to the place of use must be observed.

1. Introduction

1.4

Duty of the operating authority

The operating authority undertakes to ensure that the Robotino® is used only by persons who:

- are familiar with the basic regulations regarding operational safety and accident prevention and who have received instructions in the handling of the Robotino®,
- have read and understood the chapter on safety and the cautionary notes in this manual.

Safety-conscious working of the persons should be regularly vetted.

1.5

Duty of trainees

Prior to commencing work, all persons assigned to working on the mobile robot system Robotino® have a duty to:

- read the chapter on safety and the cautionary notes in this manual and,
- observe the basic regulations regarding operational safety and the prevention of accidents.

1.6

Risks involved in dealing with the mobile robot system

The Robotino® is designed according to state of the art technology and in compliance with recognised safety regulations. However when using the system there is nevertheless a risk of physical or fatal injury to the user or third parties or of damage being caused to the machinery or other material assets.

The Robotino® is to be used only:

- for its intended purpose and
- in an absolutely safe condition.

Faults impairing safety must be rectified immediately!

1. Introduction

1.7

Warranty and liability

In principle all our „Terms and Conditions of Sale“ apply. These are available to the operating authority upon conclusion of the contract at the latest. Warranty and liability claims for persons or material damage are excluded if these can be traced back to one or several of the following causes:

- Use of the Robotino® not in accordance with its intended purpose
- Incorrect assembly, commissioning, operation and maintenance of the Robotino®
- Operation of the Robotino® using faulty safety equipment or incorrectly fitted or non operational safety or protective devices
- Non observance of notes in the manual regarding transport, storage, assembly, commissioning, operation, maintenance and setting up of the Robotino®
- Unlawful constructional modifications on the Robotino®
- Inadequate monitoring of components subject to wear
- Incorrectly carried out repairs
- Catastrophies as a result of foreign bodies and vis major.

Festo Didactic herewith rules out any liability for damage or injury to trainees, the training company and/or other third parties which may occur during the use/operation of the system other than purely in a training situation, unless such damage has been caused intentionally or due to gross negligence by Festo Didactic.

1.8

Intended use

This system has been developed and produced exclusively for vocational and further training in the field of automation and technology. The training authority and/or the instructors is/are to ensure that trainees observe the safety precautions described in the manual provided.

The use of the system for its intended purpose also includes:

- following all advice in the manual and
- carrying out inspection and maintenance work.

2. Notes on safety

General

- Trainees must only work on the Robotino® under the supervision of an instructor.
- Observe the data in the data sheets for the individual components, in particular all notes on safety!

Electrics

- Electrical connections are to be wired up or disconnected only when power is disconnected!
- Use only low voltages of up to 24 V DC.

Mechanics

- Move the Robotino® only by carefully seizing it by its handles.
- Securely mount all components on the chassis.
- No manual intervention unless the Robotino® is at rest.

3. Technical data

Parameter	Value
Voltage supply	24 V DC, 4.5 A
Digital inputs	8
Digital outputs	8
Analogue inputs	8 (0 – 10 V)
Relais outputs	2

4. Transport/Unpacking/Scope of delivery

Transport

The mobile robot system Robotino® is delivered in two Systainers.

The Systainers must be secured against tipping over or dropping.


The carrier and Festo Didactic are to be notified immediately of any damage caused during transport.

Unpacking

Carefully remove the padding material in the Systainer when unpacking the Robotino®. When unpacking the Robotino®, make sure that none of the Robotino® assemblies have been damaged.

Note

Always seize the Robotino® by the two handles to avoid any damage to the command bridge, electronics and collision protection sensor.

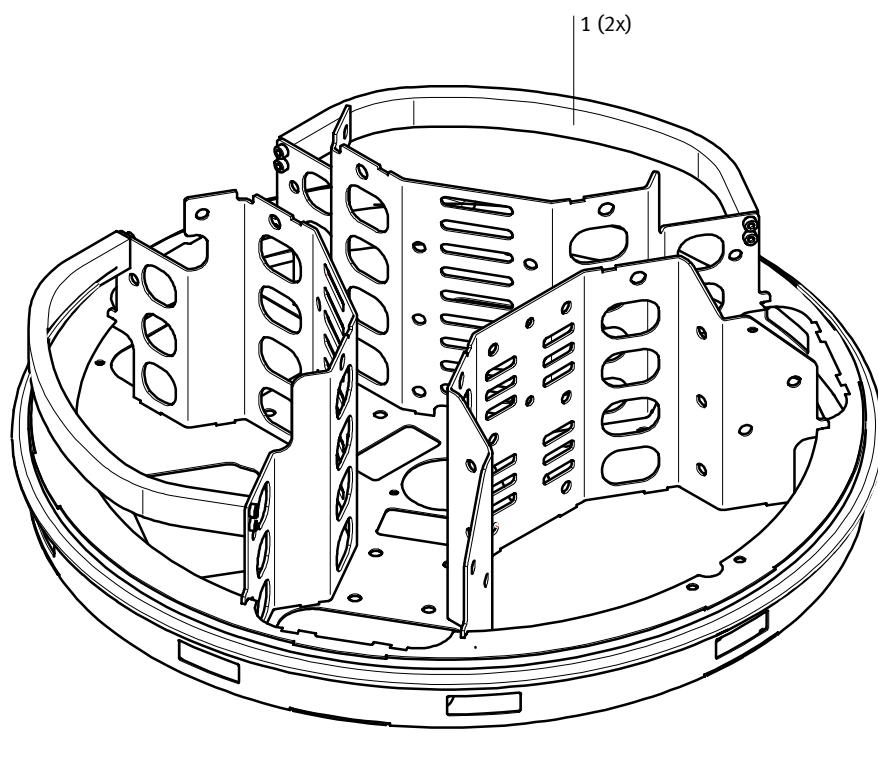Handle (1)

Check the Robotino® for any possible damaged once unpacked. The carrier and Festo Didactic are to be notified immediately of any damage.

Scope of delivery

Check the scope of delivery against the delivery note and the order. Festo Didactic must be notified immediately of any discrepancies.

5. Design and function

5.1

Robotino®

Robotino® is a fully functional, high quality mobile robot system with omnidirectional drive. The three drive units allow for motion in all directions – forward, backward and sideways – and the robot can be turned on the spot as well. It is also equipped with a webcam and several types of sensors, analogous to distance measurement, for example binary for collision protection and digital to check the actual speed. This assures that all of the wide ranging demands placed upon systems of this type are fulfilled.

The system can be placed into service immediately – without a PC.

Robotinos® controller consists of an embedded PC with a compact flash card, to which several demo applications and the operating system (Linux) have been installed. The demo applications can be started directly at the Robotino® control panel.

Robotino® can be programmed with Robotino®View software at a PC via wireless LAN. Robotino®View is capable of transmitting signals to the motor controller, as well as displaying, changing and evaluating sensor values. Robotino® can even be programmed during actual operation with Robotino®View.

Linux and C++ APIs are also available for programming Robotino®.

The webcam makes it possible to display and evaluate a live camera image with the help of Robotino®View. In this way, applications such as path and object tracking can be implemented.

Robotino® is autonomous. Numerous sensors, a camera and a high performance controller provide the system with the necessary “intelligence”.

The controller can be directly accessed via wireless LAN (WLAN). When correctly programmed, Robotino® completes assigned tasks autonomously.


Additional sensors and actuators can be connected via an I/O interface.

5. Design and function

5.2


Chassis and command bridge

The chassis consists of a laser welded, stainless steel platform.

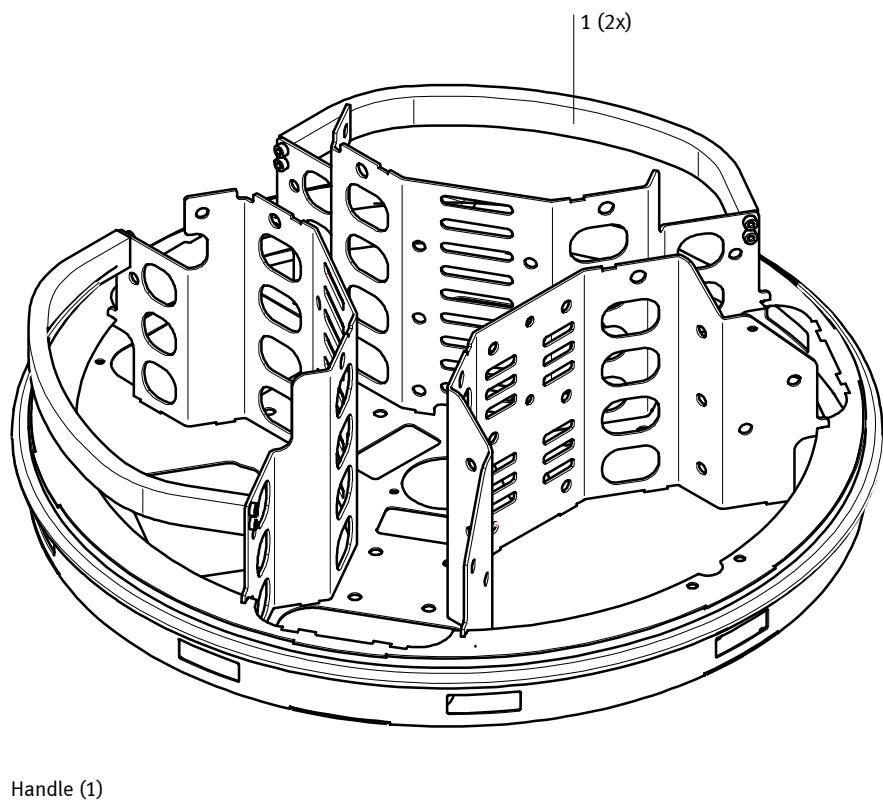
5. Design and function

The rechargeable batteries, the drive units and the camera are mounted to the chassis, where the distance measuring sensors and the anti-collision sensor are also located. The chassis offers additional space and mounting options for other attachments, sensors and/or actuators.

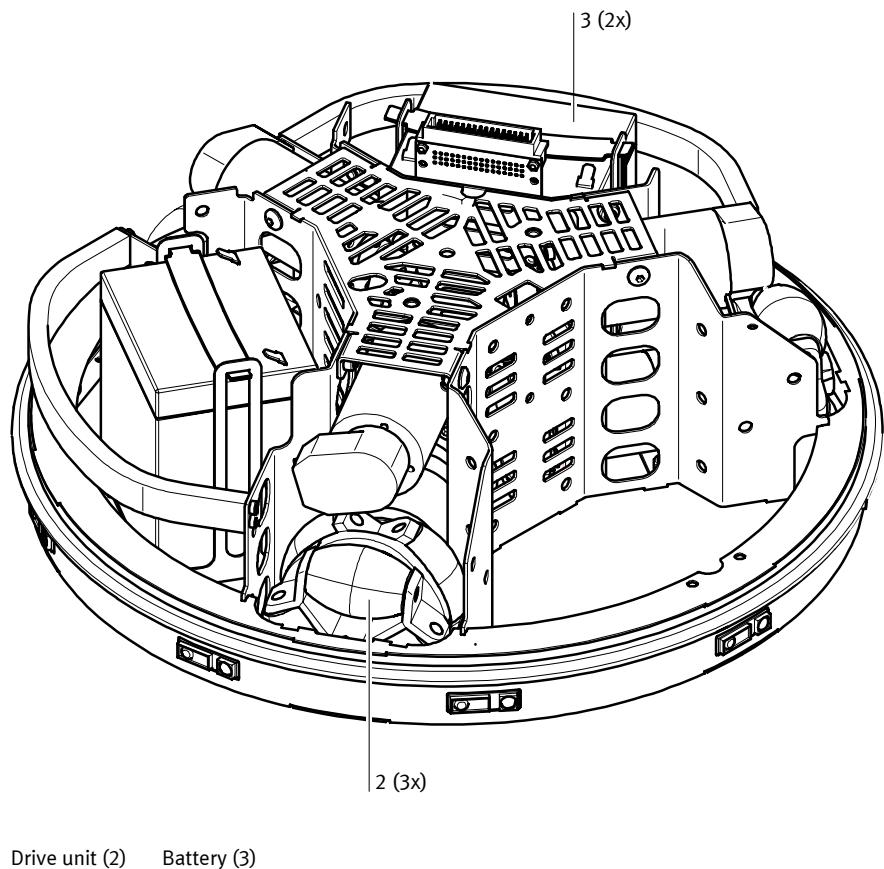5. Design and function

Highly sensitive system components such as the controller, the I/O module and the interfaces are situated on the command bridge. The command bridge is connected to the other system modules by means of a plug connector.

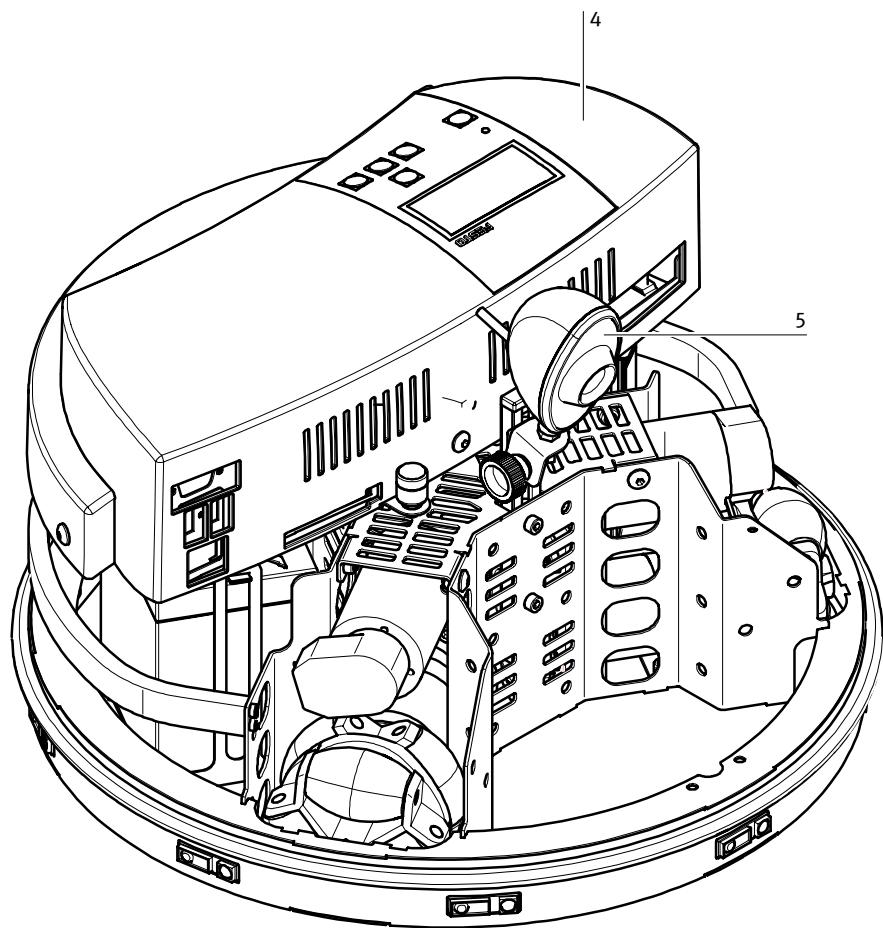Command bridge (4) Camera (5)


5.3

The drive unit module

Robotino® is driven by 3 independent, omnidirectional drive units. They are mounted at an angle of 120° to each other.

Each of the 3 drive units consists of the following components:

- DC motor
- Gear unit with a gear ratio of 16:1
- All-way roller
- Toothed belt
- Incremental encoder

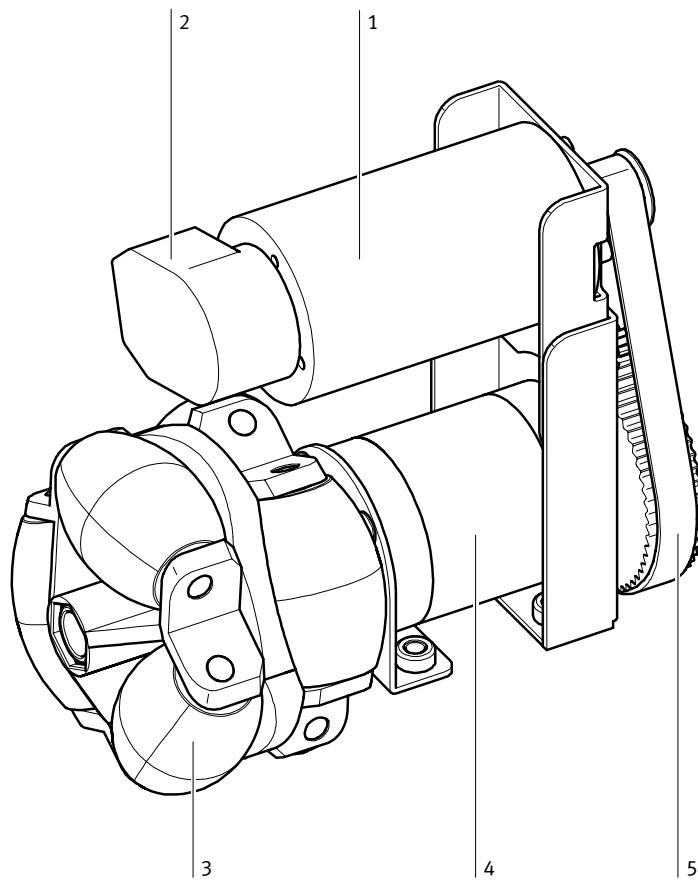Motor (1)

toothed belt (5)

incremental encoder (2)

all-way roller (3)

gear unit (4)

All of the individual components are attached to the mounting flange at the rear. Together with the flange at the front, the drive unit is secured to the chassis with screws. Correct positioning of the drive units in relationship to each other is thus assured.

5. Design and function

Actual motor speed can be compared with desired speed by means of the incremental encoder, and can then be regulated with a PID controller via the I/O circuit board.

Motor performance data

DC motor (GR 42x25)	Unit of measure	
Nominal voltage	V DC	24
Nominal speed	RPM	3600
Nominal torque	Ncm	3.8
Nominal current	A	0.9
Starting torque	Ncm	20
Starting current	A	4
No-load speed	RPM	4200
No-load current	A	0.17
Demagnetisation current	A	6.5
Mass moment of inertia	gcm ²	71
Motor weight	gr.	390

Gear unit

Planetary gear unit (PLG 42 S)	
Single-stage, Nm:	3.5
Single-stage, i:	4 :1 – 8 :1
2-stage, Nm:	6
2-stage, i:	16 :1 – 64 :1
3-stage, Nm:	14
3-stage, i:	100 :1 – 512 :1

All-way roller

All-way roller, driven (ARG 80)	
Diameter Ø	80 mm
Maximum carrying capacity	40 kg

The all-way roller is set into motion in a given direction by means of its drive axle, and is also capable of travelling in any desired direction when powered by rollers. As a result of interaction with the other two drive units, travel is thus possible in a direction which deviates from the direction of the respective drive.

5.4 The camera module

Robotino® is equipped with a camera system. Its height and inclination can be adjusted. The camera makes it possible to display live images with the help of Robotino®View. Robotino®View also offers a number of image processing options, which can be used to evaluate camera images for the Robotino® controller. A segmenter locates surfaces of like colour within a given image, and can determine the position and size of any segment. Lines in video images can be detected. Results can be utilised for pinpointing objects, as well as for path and object tracking.

Technical Specifications	
Image Sensor	Colour VGA CMOS
Colour Depth	24 Bit True Colour
PC-connection	USB 1.1
Video resolutions	160 x 120, 30fps (SQCGA) 176 x 144, 30fps (QCIF) 320 x 240, 30fps (QVGA) 352 x 288, 30fps (CIF) 640 x 480, 15fps (VGA)
Still Image Resolutions	160 x 120 (SQCGA) 176 x 144 (QCIF) 320 x 240 (QVGA) 352 x 288 (CIF) 640 x 480 (VGA) 1024 x 768 (SVGA)
Still Capture Format	BMP, JPG

Note

If at all possible, plug the camera cable into the USB port on the right-hand side, because this reduces the risk of destruction of the cable or the camera which results from the cable protruding beyond the housing.

5.5

The controller unit


The controller unit can be set up flexibly through the use of several plug-in modules.

The Robotino® controller unit consists of 3 components:

- PC 104 processor, compatible with MOPSlcdVE, 300 MHz, and Linux operating system with real-time kernel, SDRAM 128 MB
- Compact flash card (256 MB) with C++ API for controlling Robotino®
- Wireless LAN access point

The controller unit is equipped with the following interfaces:

Ethernet, 2 ea. USB and VGA. These are used to connect a keyboard, a mouse and a monitor screen. The operating system and the C++ library can thus be accessed without a PC, if a WLAN connection is not possible or not desired. The Ethernet connection cannot be used with the basic version.

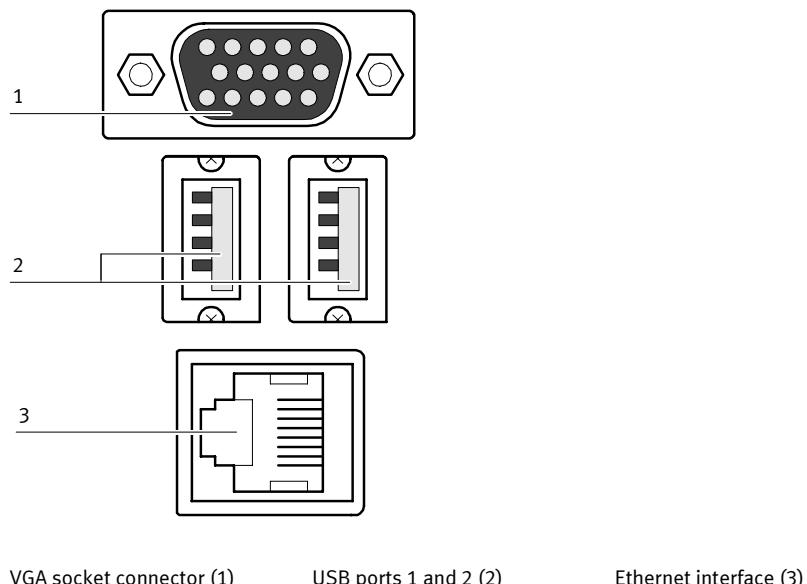5.6

I/O circuit board module

The I/O circuit board establishes communication between the controller unit and the sensors, the drive units and the I/O interface included with Robotino®.

Each of the motors in the individual drive units are controlled by a PID controller. Each motor can be controlled individually.

Signals from the step encoder, all permanently installed sensors and all sensors and actuators which are connected to the I/O interface are forwarded to the controller unit or the additional actuators.

5.7

Power supply / battery charger

Power is supplied by two rechargeable 12 V batteries with a rating of 4 Ah. Both of the rechargeable batteries are mounted to the chassis. Robotino® is supplied with 2 additional batteries and a battery charger. Two batteries can thus be recharged, while the other two are in use.

5.8

Sensors

Sensors for measuring the distance to objects and detecting motor speed have been integrated into Robotino®. An anti-collision sensor mounted around the circumference of the chassis signals contact with objects.

Infrared distance measuring sensors

Robotino® is equipped with 9 infrared distance measuring sensors which are mounted in the chassis at an angle of 40° to one another. Robotino® can scrutinise all surrounding areas for objects with these sensors. Each of the sensors can be queried individually via the I/O circuit board. Obstacles can thus be avoided, clearances can be maintained and bearings can be taken on a selected target. The sensors are capable of accurate or relative distance measurements to objects at distances of 4 to 30 cm. Sensor connection is especially simple including just one analogue output signal and supply power. The sensors' evaluation electronics determine distance and read it out as an analogue signal.

Incremental encoder

The actual speed of each individual motor is measured in RPM by the incremental encoder. If actual motor speed deviates from the setpoint, it can be adjusted to match the target value by means of a PID controller whose parameters are configured with the help of Robotino® View software.

The anti-collision sensor

The anti-collision sensor is comprised of a switching strip which is secured around the entire circumference of the chassis. A switching chamber is located inside a plastic profile. 2 conductive surfaces are situated inside the chamber, between which a given clearance is maintained. These surfaces are short circuited when even minimal pressure is applied to the strip. A reliably recognisable signal is thus transmitted to the controller unit. Collisions with objects at any point on the housing are detected and, for example, Robotino® is brought to a standstill.

5. Design and function

Analogue inductive proximity sensor

The inductive proximity sensor is supplied as an additional component. It serves to detect metallic objects on the floor and is used for continuous-path control. It reads out signals of varying strength depending upon whether it is located in the middle or at the edge of the metal strip. Path tracking can thus be controlled in a differentiated fashion.

The inductive proximity sensor must be attached to the mounting furnished for this purpose, and must be connected to the I/O interface.


Technical data	
Operating voltage	15 – 30 V DC
Output voltage	0 – 10 V
Type	SIEA-M12B-UI-S
Part no.	538292
Diameter	M12
Sensing range	0 to 6 mm
Mounting	Quasi embeddable
Switching frequency	1000 Hz
Ambient temperature	-25 – +70° C
Protection	IP 67
Housing material	Chrome-plated brass
Max. tightening torque	10 Nm
Reproducibility	0.01 mm

5. Design and function

Diffuse sensors

Path tracking can also be implemented with the two included diffuse sensors. Flexible fibre-optic cables are connected to a fibre-optics unit which works with visible red light. Reflected light is detected. Different surfaces and colours produce different degrees of reflection. However, gradual differences in reflected light cannot be detected.


The sensors must be attached to the mountings furnished for this purpose, and must be connected to the I/O interface.

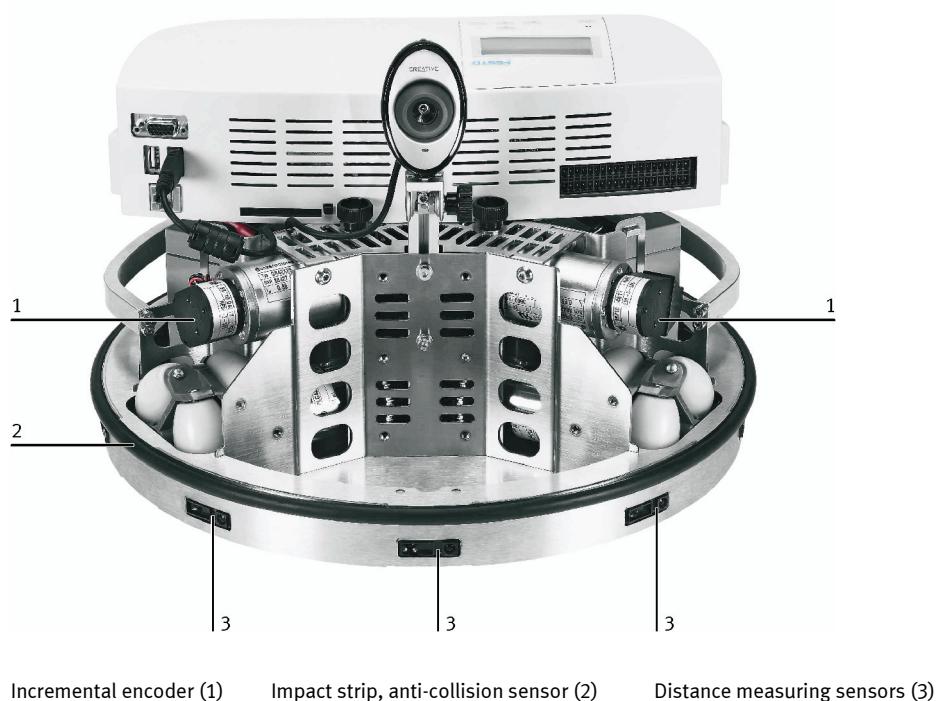Sensor positions and their designations in Robotino® View

The distance measuring sensors and the incremental encoder must be uniquely identified in order to be able to address them with Robotino® View in a targeted fashion. The respective designations are included in the figure below. IR1 is addressed as “Distance 1” in the software. In a like fashion, the other infrared sensors are designated “Distance 2” through “Distance 9”.

The incremental encoders are assigned to their respective drive units.

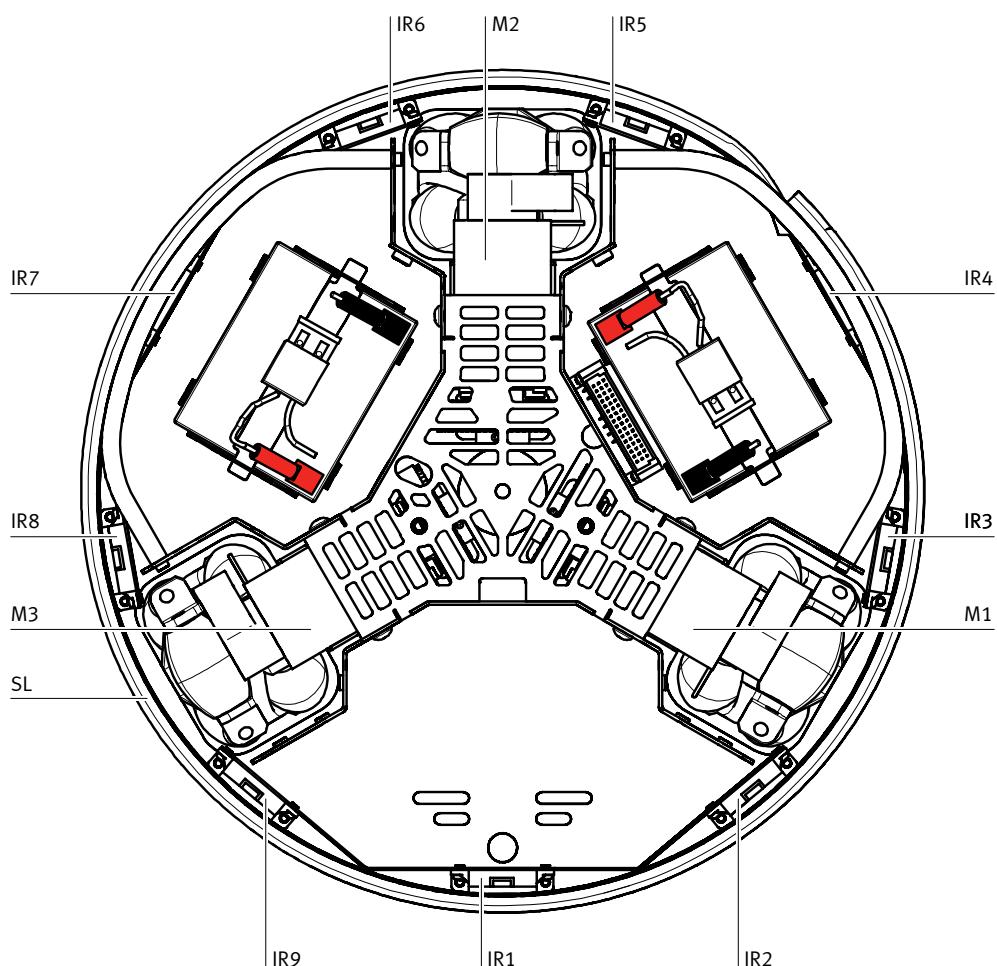Sensor assignments, IR1 to IR9: distance measuring sensors (1 – 9)

M1 to M3: motors (1 – 3)

SL = impact strip, anti-collision sensor


5.9

Membrane keypad and display

A membrane keypad and a display are located on the housing top, by means of which various options can be selected, information can be queried, and included programs can be started.

Example of an initial display

Display text	Description
Robotino®	
172.26.1.1	Robotino® IP address
V1.0	Software version

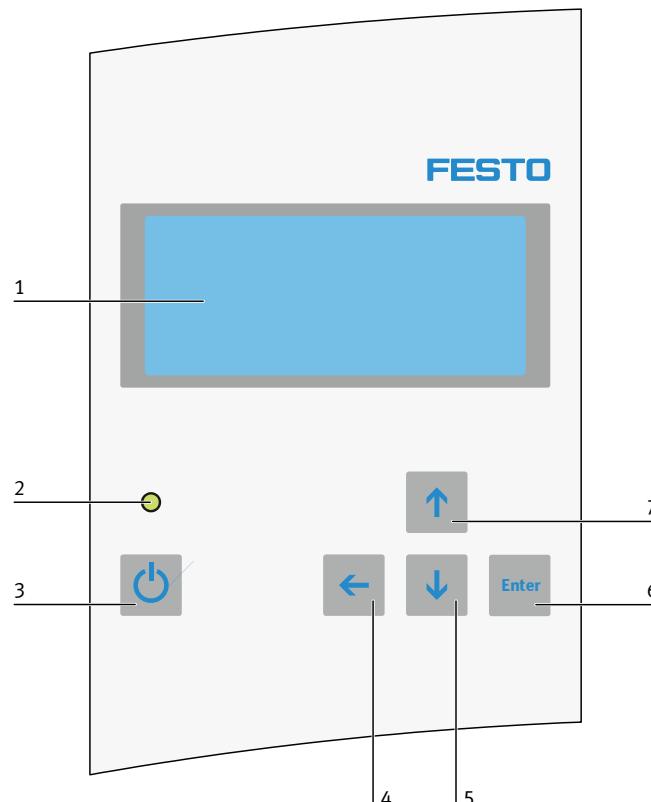Display (1) LED (2) on/off (3) shift up one menu level (4)
 scroll down one selection (5) acknowledge selection (6) scroll up one selection (7)

5.10
Wireless LAN access point

The wireless LAN access point is a component which allows for communication with the robot via a network address.

- The access point features low current consumption. Power supply via the USB port is possible.
- The access point complies with the following standards: IEEE 802.11g and 802.11b.
- It allows for transmission speeds of up to 54 Mb per second for 802.11g, and 11 Mb per second for 802.11b with a large transmission range (up to 100 m inside of buildings).
- It assures reliable network security with WEP encryption and WPA-PSK function.
- It's fast and simple to configure via the web management utility.

5.11
The compact flash card

The controller unit is equipped with a slot into which a PC card has already been inserted. This PC card contains the operating system, the functions libraries and the included programs. Updates can be easily installed by simply replacing the PC card. The slot for the PC card is located to the right of the controller unit interfaces.


5.12

The I/O interface

The I/O interface makes it possible to connect additional sensors and actuators.

They are connected by means of an included plug.

- 8 analogue inputs (0 to 10 V) (AIN0 to AIN7)
- 8 digital inputs (DI0 to DI7)
- 8 digital outputs (DO0 to DO7)
- 2 relays for additional actuators (REL0 and REL1). The relays can be connected as NC, NO or of CO contacts.

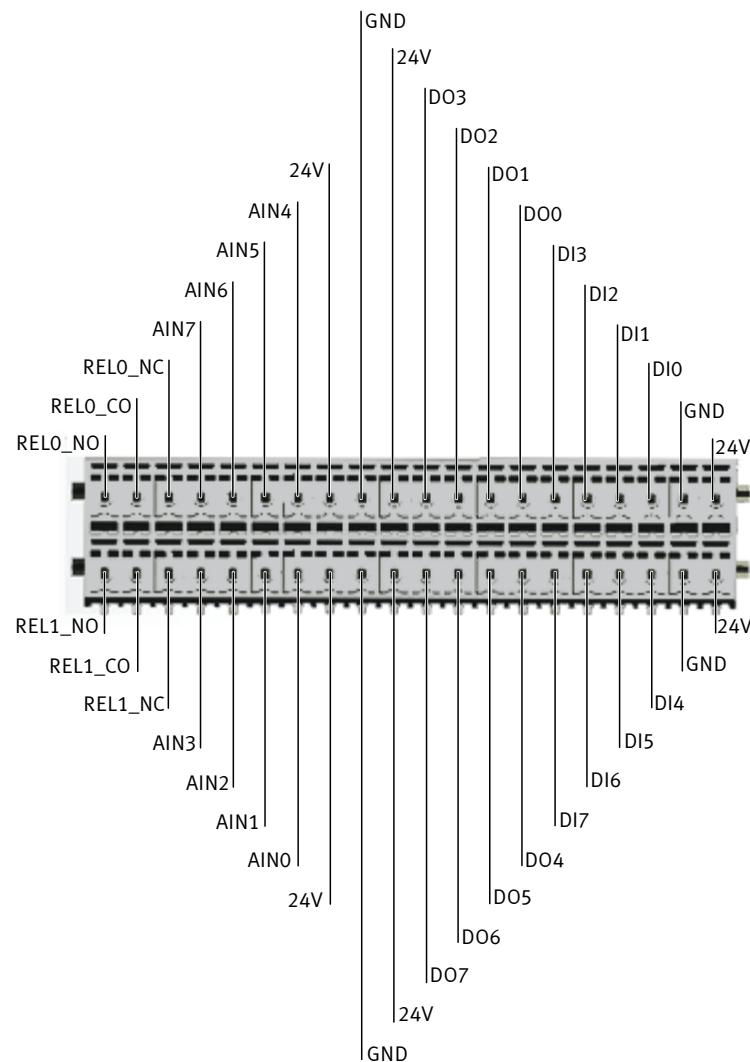I/O interface pin allocations

6. Commissioning

6.1

Initial steps

Always jack up the Robotino® if you do not want to work with the Robotino® or would like to test individual functions! You can then easily check the functioning of sensors or activate actuators without fear of a collision. In this way, the Robotino® is always close at hand and can remain connected to the charging station.

With the supplied device, you can jack up the Robotino® such that the Robotino® wheels are freely movable and do not come into contact with the ground. To do so, place the jacking device onto a stable and level surface with the rubber pads pointing upwards. Then position the Robotino® onto the rubber pads. The Robotino® must be positioned onto the jack such as to ensure that the drive pinions are not located on any of the rubber pads.

Important notes

Always hold the Robotino® at the chassis when transporting it. Lifting the Robotino® from the command bridge can cause extensive damage to the sensitive electronics. The front of the Robotino® is where the distance sensor 1 and the mounting facility for the inductive sensor are located.

Any obstacles must extend as far as the ground since otherwise the collision protection sensor will not respond. The same applies in the case of the distance sensors. This is particularly important since otherwise damage may be caused to the Robotino®.

Next, connect the battery charger to the appropriate cable of the batteries. The cable is located behind the drive unit M1. The plug is under the ring of the base group. Carefully pull this out and connect it to the charger by inserting the plug into the socket, making sure that the latching clamps of the charger cable is clipped in. To release the cable connection, you need to press the latching clamp and pull out the plug at the same time.

6.2

Switching On and Off

Now switch on the Robotino®.

Press the On/Off button until the LED lights up.

The display switches on.

On the display, 2 bars are displayed across the entire display width.

The PC of the Robotino® is now booting up.

After approx. 30 seconds, the start display appears on the display.

The display shows for example the following: Robotino®
172.26.1.1

The last line shows a bar indicating the state of charge of the batteries and displaying Robotino® version: V 1.0

The Robotino® is now ready for operation.

If none of the keys are actuated for 10 seconds, the illumination of the display is switched off in order to keep current consumption as low as possible during operation. To re-activate the display, press one of the arrow keys to switch on again.

Note Do not press the Enter key for this so as to avoid an unwanted start of, for example, a demo program!

Switching off

Press the On/Off button until the LED is extinguished and then release the button. The Robotino® is not switched off until the button is released.

6.3

The display functions

Main menu

Press the Enter key to access the main menu. In its original status the display is in English.

It contains the following menu items:

- Languages
- State of charge
- DEMOs
- Network

Using the Up and Down arrow keys, you can move the left (>) sentinel up or down. Press the Enter key to call up the relevant menu item or start the tagged demo program.

Press the Left arrow key to exit from a menu to the next higher-order menu.

Languages

deutsch

english

français

español

State of charge

This display informs you of the actual state of charge of the accumulators, the actual voltage and current.

Bar display for state of charge

Voltage: 23 V

Current: 1.0 A

DEMOs

Circle

Forward

Quadrangle

Roaming

Follow line

Network

This menu displays the current IP address which can be changed. The cursor (head of the arrow) is moved to the right using the Enter key and to the left using the Left arrow key. The value of the individual numbers can be increased or reduced using the Up/Down arrow keys. To exit the menu again, move the Left arrow key to the first digit. To return to the main menu again, press the Left arrow key. Press the Enter key again at the last digit of the IP address to access the next display, the subnet mask. Proceed as in the case of the IP address. When you are on the last digit of the subnet mask your changes to the IP address and subnet mask will be copied by pressing the Enter key again.

6.4

Testing the demo programs

First, check the functioning of the demo programs in the jacked up condition in order to familiarise yourself with their functions. Select a demo program and start this by pressing the Enter key.

Stop the programs by briefly pressing any one of the keys on the touch-sensitive keyboard or by pressing the bumper. The Robotino® immediately stops in either case.

The Circle and Forward programs are subject to a time limit in that they automatically stop the Robotino® after 10 seconds.

Circle

The Robotino® moves along a circular path
Required space: A rectangle of a 1m x 1m

Time limit: 10 seconds
End program: See above

Forward

The Robotino® moves straight ahead
Required space: A straight distance of 1 m
Time limit: 10 seconds
End program: See above

Quadrangle

The Robotino® moves along a rectangle and rotates around its own axis
1 m x 1m
Time limit: None
End program: See above

Roaming

The Robotino® moves straight ahead and evades any obstacles in the way.
Depending on the position of the obstacles, the Robotino® changes its direction of travel in different ways.
Required space: Any
Time limit: None
End program: See above.

Comment

Distance sensors 1, 2 and 9 are used to detect obstacles. This is why the program cannot be stopped via the bumper within the coverage range of these sensors as the Robotino® takes evasive action.

Follow line

The camera of the Robotino must be connected to the control unit via the USB interface for this program. Draw a line which is at least 5 cm wide on the floor, for example using the adhesive tape. The line must be in red, black or blue. Set the Robotino® onto the line so that it can detect the line and then start the program. The Robotino® will follow the marked line.

6. Commissioning

6.5

Executing demo programs on solid ground

The ground on which the Robotino® is to travel must be horizontal and level. In this way, the desired movements can be reliably executed. The ground must be clean and dry in order to avoid any damage to the mechanical and electronic components. You must release the cable of the charger prior to starting the program and stow it away in the chassis.

Notes

Depending on the condition of the ground, the required space specified above may vary. You must therefore select a sufficiently large clear area. The required space specified above represents a minimum requirement.

6.6

Controlling the Robotino® with Robotino®View

A WLAN connection must be established to enable you to control the Robotino® with Robotino®View.

6.7

Setting up a WLAN connection

Since different WLAN hardware and software is available, it is not possible to describe the correct set-up and preparation of your WLANs in detail.

Activate your WLAN. Your WLAN setting must permit the assignment of an IP address to the WLAN. Only then can the access point in the Robotino® establish a connection to your network.

As a precaution, place the Robotino® onto the jacking device. Then switch on the Robotino® and wait for the booting up process, which you can view on the display. Make a note of the IP address displayed. To be on the safe side, a label with the IP address is affixed to the underside of the command bridge.

Let your WLAN search for available networks. A network with the name Robotino® x.x now appears in the list of available networks. If necessary, establish a connection to this network if this has not been effected via your network software.

The following network settings are a prerequisite for a successful connection:

Automatically assign network code (SSID)

Automatically source IP address

Both settings must be switched on to be able to establish a connection to the Robotino®.

6.8

Testing the WLAN connection

You can test the WLAN connection using 2 simple DOS commands or the Robotino®View software.

The ping IP address command

Start the MS-DOS prompt, which is in the start menu under Programs/Accessories. Click onto "Prompt".

Alternatively, you can enter the command "cmd" under "execute" in the start menu. Now enter the „ping“ command in the prompt, together with the IP address you have read off the Robotino® display, e.g. "**ping 172.26.1.1**"

You will receive the following messages if a WLAN connection exists:

```
C:\>ping 172.26.1.1
```

Ping is executed for 172.26.1.1 with 32 bytes of data:

```
Response from 172.26.1.1: Bytes=32 time=4ms TTL=64
Response from 172.26.1.1: Bytes=32 time=2ms TTL=64
Response from 172.26.1.1: Bytes=32 time=3ms TTL=64
Response from 172.26.1.1: Bytes=32 time=6ms TTL=64
```

Ping statistics for 172.26.1.1:

Packages: Sent = 4, received = 4, lost = 0 (0% loss),

Approx. time specifications in milliseconds:

Minimum = 2ms, maximum = 6ms, mean value = 3ms

The following message will be displayed in the absence of a connection to the Robotino®:

Ping is executed for 172.26.1.1 with 32 bytes of data:

Timeout of request.

Timeout of request.

Timeout of request.

Timeout of request.

Ping statistics for 172.26.1.1:

Packages sent = 4, received = 0, lost = 4 (100% loss),

The command ipconfig /all

This command enables you to check the settings and statuses of all your network connections. Start the MS-DOS input request, which is available in the start menu under programs/accessories. Click onto "input request".

Alternatively, you can enter the command „**cmd**“ in the Start menu under Execute.

Enter the command "**ipconfig /all**" in the input request ">" and press the Enter key.

Once the connection to the Robotino® is established, you will receive the following information. However, different options are possible here depending on the type of your WLAN adapter. Find the setting for wireless network connection under the various network settings.

Ethernet adapter wireless network connection:

Connection-specific DNS-Suffix:

Description.....: Dell Wireless 1470 Dual Band WLAN Mini-PCI card

Physical address: 00-14-A5-44-A4-D6

DHCP activated: Yes

Autoconfiguration activated....: Yes

IP address: 172.26.201.1

Subnetmask.....: 255.255.0.0

Standard gateway.....:

DHCP-Server: 172.26.101.8

Lease received: Friday, 17. February 2006 13:51:14

Lease expires: Friday, 24. February 2006 13:51:14


6. Commissioning

Checking the WLAN connection with Robotino® View:

Start the Robotino® View software.

Enter the IP address of the Robotino® in line (1). Click onto the aerial symbol (2) to the right of the input line

Once a connection is established, the message “connected to 172.26.1.1.” is displayed in the status bar (3) on the bottom left.

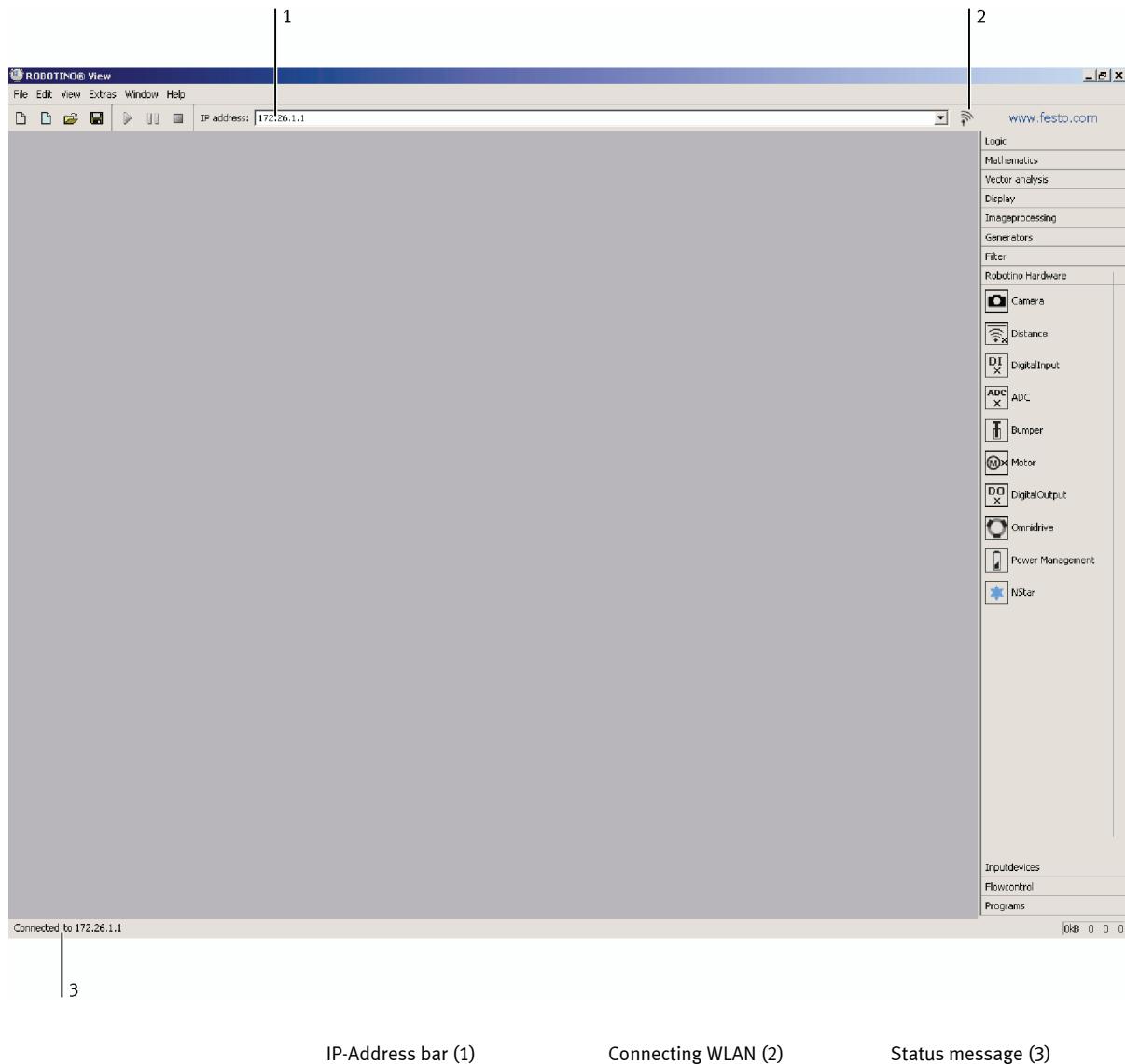6. Commissioning


If there is no connection to the Robotino®, a message is output after a certain time (45 s) and the message “connection refused” is displayed in the status bar on the bottom left. The aerial symbol is animated whilst the software is trying to establish a connection.

If the connection has been successfully set up, you can access all the software functions and all Robotino® units directly.

By way of a simple example, now access the webcam of the Robotino®. To do so, connect the webcam to the command bridge via the USB interface.

- Open a new worksheet by clicking onto the symbol (1) in the toolbar.
- Now open the "Robotino® hardware" folder (2) by clicking it.
- Drag the camera symbol (3a) onto the desktop by holding down the mouse button.
- Click onto the green arrow in the toolbar (4) to start the program
- Open a window (5) with the actual camera image by double clicking the camera symbol on the desktop (3b). By reducing the desktop, you can enlarge the window showing the camera image.

6. Commissioning

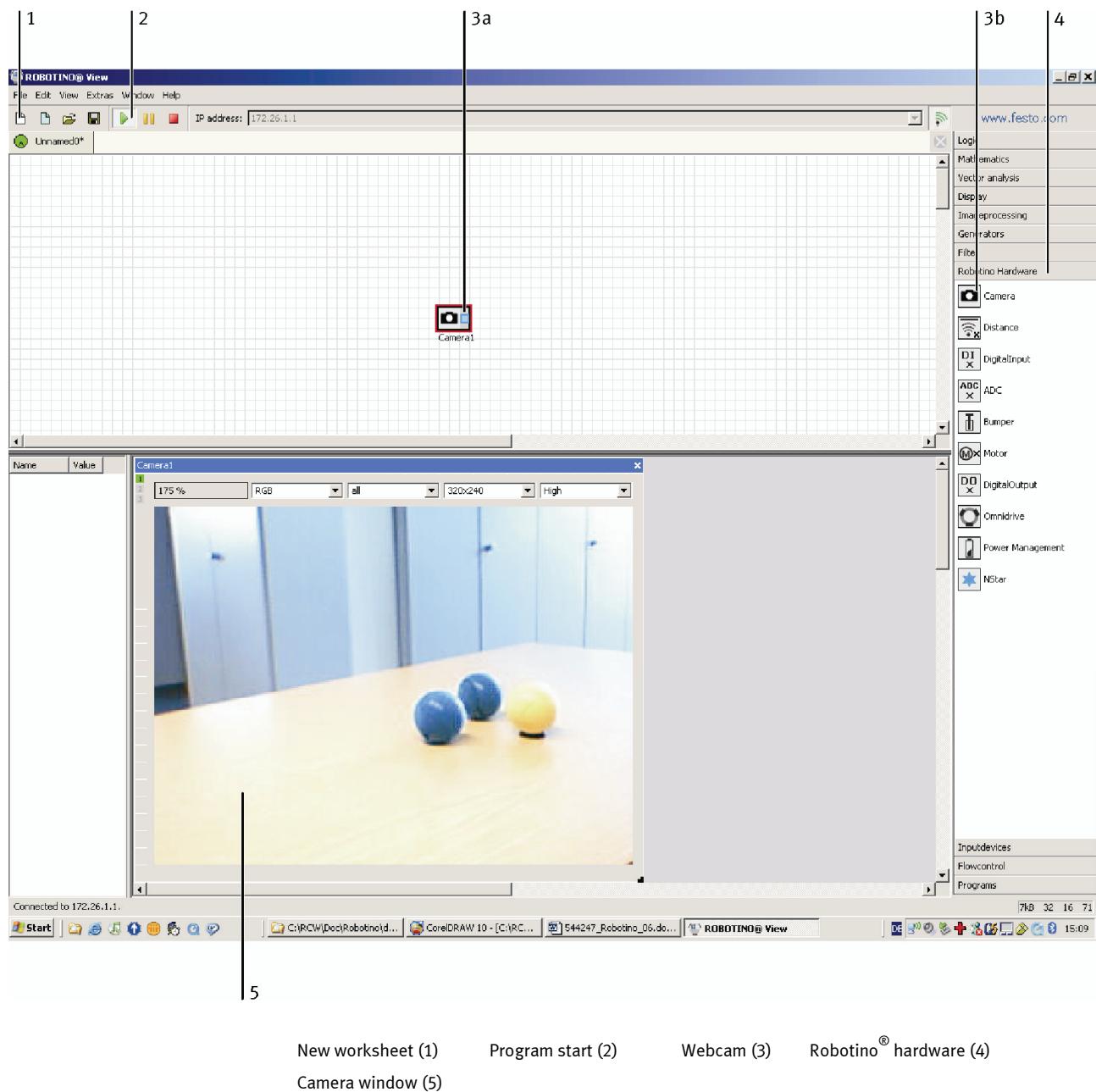
6.9

Working with several Robotinos®

Working with one Robotino®

The Robotino® has its own WLAN server. When operating one Robotino® you will therefore only need a PC which can establish a WLAN connection. In this mode of operation, the WLAN server of the Robotino® is in the AP (Access Point) mode.

Advantage: It is not necessary to effect a configuration of the AP and network respectively. All that is required is WLAN capable PC.

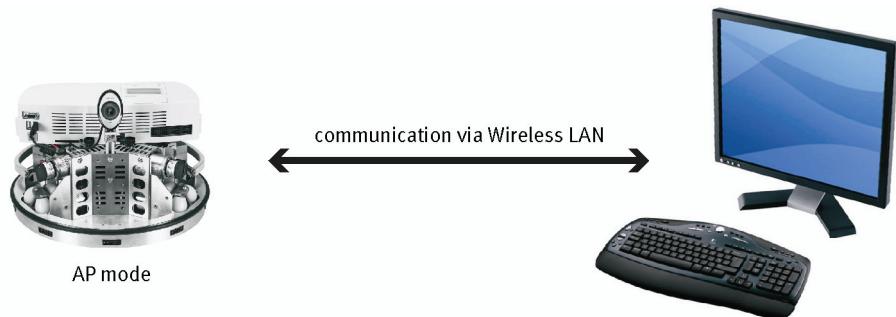Working with three to four Robotinos®

If three or four Robotinos® are to be controlled simultaneously, the WLAN server of the Robotino® can be operated in the AP mode described above. One advantage of this is that all Robotinos® can share the same IP address, since each one forms its own network. It is therefore possible to access different Robotinos® from each PC without the need to change or establish the IP address.

The disadvantage of this mode of operation is that different WLAN networks can cause collisions if their transmission channels are too close together. A maximum of 11 of these channels is available for a WLAN network. For reasons of safety it is advisable that at least three non active channels should remain free between two active channels. In this way a maximum of 4 Robotinos® can be safely operated.

Advantage: Several Robotinos® can be operated independently of one another.


Disadvantage: The transmission channel must be defined for each Robotino® and for each PC.

Working with several Robotinos®, if the PCs are connected to a training network.

For this mode of operation the access point of the Robotinos® must be set to client mode. This is effected via a switch directly on the access point of the Robotino®. In this case a central WLAN access point is required which is connected to a local Ethernet network.

Advantage: Any number of Robotinos® can operate in one network.

Disadvantage: Each Robotino® requires its own IP address, which must be entered via the touch sensitive keyboard of the Robotino®.

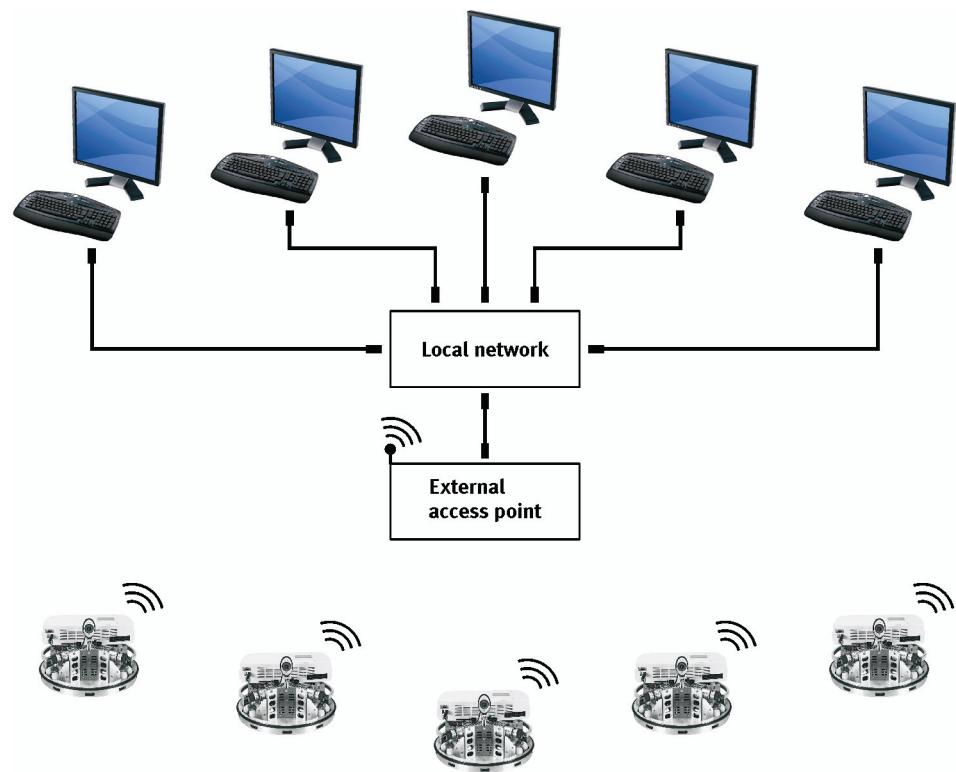Also the local network is accessible via the unencrypted access point.


Settings	Value
SSID	RobotinoAPx.1
Channel	11
Encryption	none

Working with several Robotinos® in the absence of a training network

The WLAN of the Robotinos® must be set to the AP client mode; this is effected via a switch directly on the Robotino®. A central additional WLAN server is then required.

Advantage: Any number of Robotinos® can operate within a network.

Disadvantage: Each Robotino® requires its own IP address; this must be entered via the touch sensitive keyboard of the Robotino®.

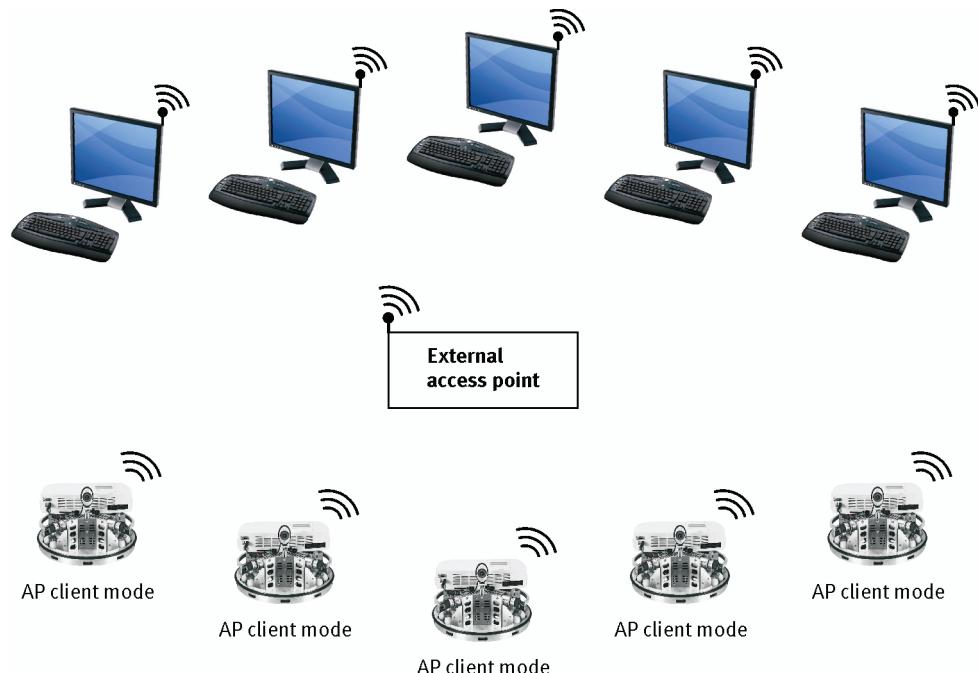Also the local network is accessible via the unencrypted access point.

6.10 Charging the batteries


Two Additional batteries are supplied with Robotino®. Batteries which have been depleted during use can be replaced and recharged while the second set is being used. Robotino® can thus be operated without interruption. In order to charge the batteries, set them next to each other such that the positive terminal of one battery (identified in red) and the negative terminal of the other battery (identified in black) are directly next to each other. Connect these two terminals with the included blue connector cable. Push the cable lugs on the ends of the cable over the battery terminals to this end. Then connect the clamps on the ends of the battery charging cable to the other two battery terminals. Connect the red clamp to the vacant positive terminal (identified in red) and the black clamp to the vacant negative terminal (identified in black). Connect the battery charging cable to the cable from the battery charger. Insert the plug on the end of the cable from the battery charger into the socket on the end of the battery charging cable to this end. Make sure that the latch on the battery charging cable snaps into place.

6. Commissioning

Notes

Always charge both batteries simultaneously. If an individual battery is charged, this will result in it being damaged.

If the rechargeable batteries are completely empty, these must be recharged again to ensure the correct functioning of the Robotino®. If the batteries are not sufficiently charged, the functioning of the drive units may be erratic.

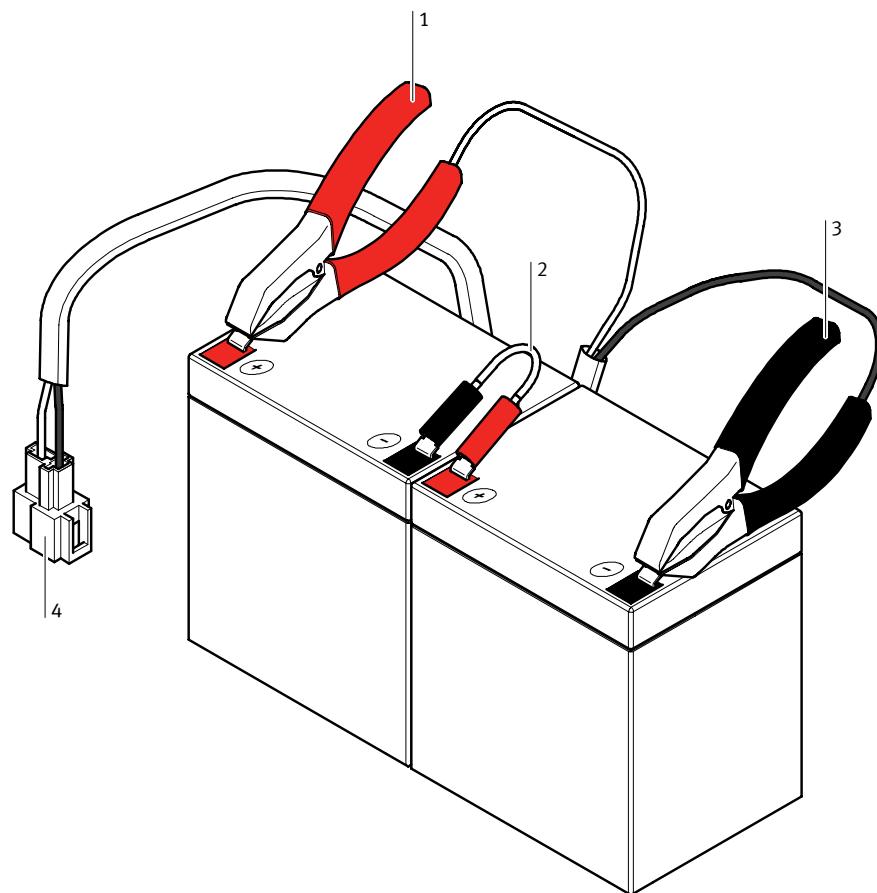Positive terminal (red) (1)
latch (4)

Blue connector cable (2)

Negative terminal (black) (3)

The latch in the plug on the end of the battery charging cable must be pressed while simultaneously pulling out the plug in order to disconnect the cables.
LED display on the battery charger:


LED colour	State
Red	power on
Yellow	charging
Green	charging completed

6.11 Replacing the batteries

In order to replace the batteries, the command bridge must first be removed from the chassis. Proceed as follows:

Required tool: 4 mm hex key

- Undo all of the plug connections on the command bridge.
- Remove the two hexagon socket head screws from the front, as well as the two hexagon socket head screws from the centring plates underneath the command bridge in order to dismantle it from the chassis.
- Pull the command bridge up and out, and set it carefully aside.
- Disconnect the cable lugs from the negative and positive battery terminals.
- Open the battery holder.
Press down both ends of the clasp (1) on the top of the battery holder to this end, and bend out the two retaining brackets (2). Set the clasp aside.
- Remove the battery (3). Tip the battery towards the outside to this end, and lift it out.

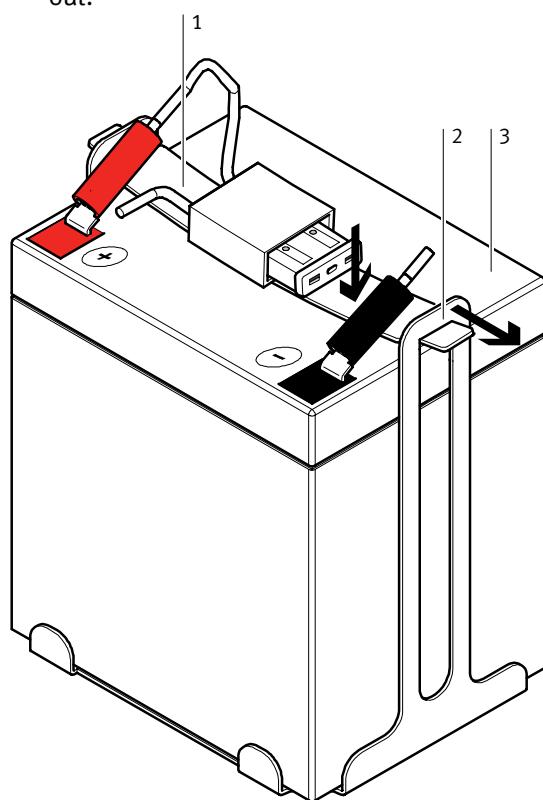Clasp (1)

Retaining brackets (2)

Battery (3)

6. Commissioning

Repeat the same procedure for the second battery.

Repeat the procedure in reverse in order to install the replacement batteries.

Observe the following points:


- Make sure than none of the battery cables are trapped underneath the batteries, or between the batteries and the battery holders.
- Both of the battery terminals must point towards the centre of the Robotino[®].

The cables must be connected to the battery terminals as follows:

- The black cable must be connected to the negative terminal on the left (black).
- The red cable which is secured to the clasp with the locking device must be connected to the positive terminal at the right (red).

Securing the command bridge

From the rear, guide the centring plates (2) on the command bridge up to the mounting brackets (3) on the chassis. Gently press the command bridge down in order to insert it into the connector plug. Mount the command bridge securely with the previously removed screws.

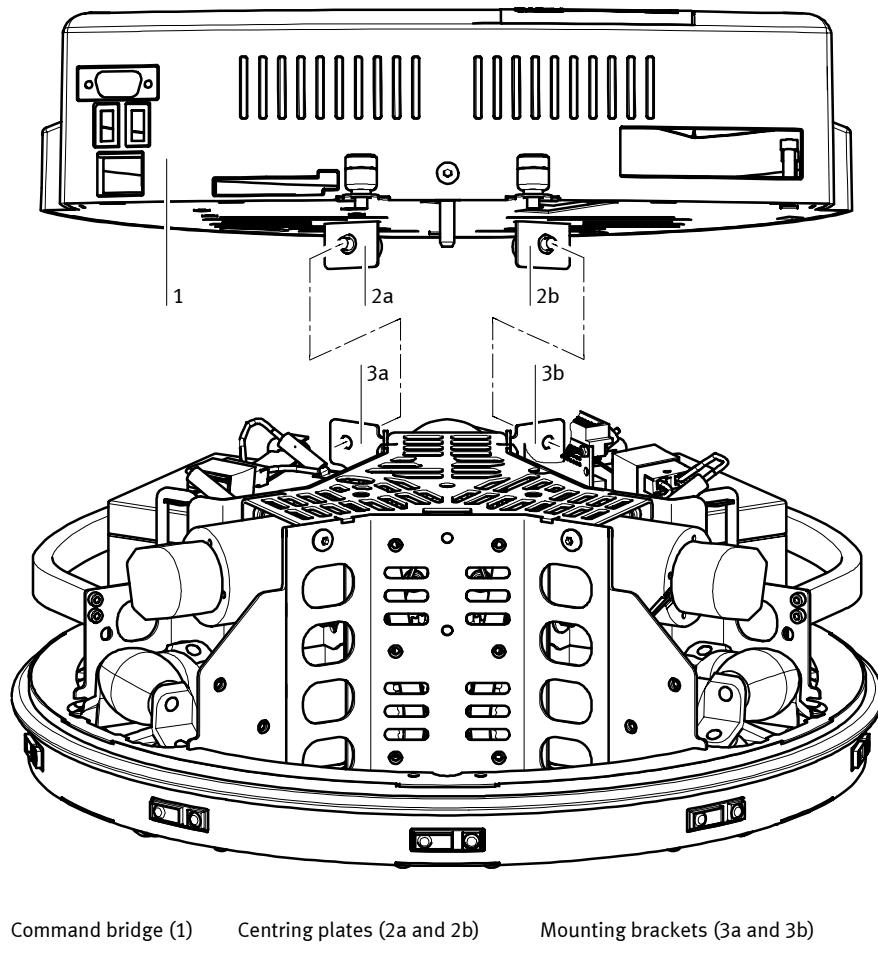6.12 Installing additional sensors

Diffuse sensor

Required tools: 2, 5 and 3 mm hex keys, 10 mm open-end wrench, 3 mm Phillips head screwdriver


All of the components for the diffuse sensor must be pre-assembled before they are mounted to Robotino® and connected.

The following components are required for each sensor:

- 1 fibre-optic unit
- 1 fibre-optic cable
- 1 retainer for the fibre-optic cable

First of all, cut the two fibre-optic cables to the required length. It is absolutely essential to use a fibre-optic cutting device to this end, in order to assure that the fibre-optic cables are no destroyed, and to ensure error-free functioning. Screw the fibre-optic head into the mounting bracket until it protrudes from the other side by approximately 2 mm. Secure it with the included nut. When the mounting brackets are attached to the chassis during a later assembly step, it must be assured that the sensors are located at the inside ends of the slots. Otherwise, the distance between the sensors cannot be changed.

Insert the free end of the fibre-optic cable into the black retainer in the fibre-optic unit. Push it all the way in until resistance is plainly perceived. Tighten the retainer with the Phillips head screw.

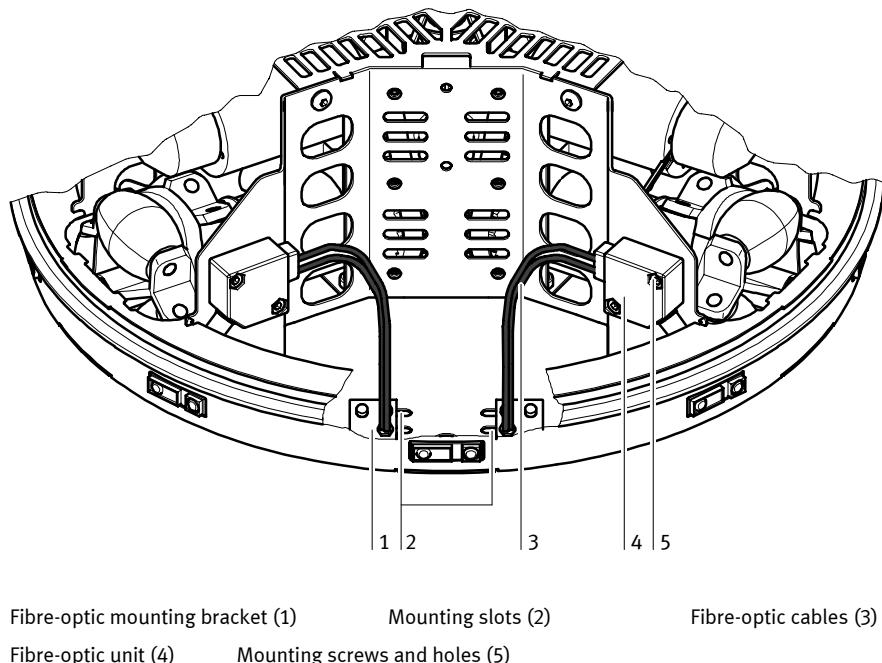6. Commissioning

Now mount the sensors to Robotinos® chassis. First of all, secure the fibre-optic mounting bracket (1) to the bottom of the chassis by inserting two screws into the slot from underneath. The slots (2) at the front of the chassis for the mounting bracket have different functions. The long slot is used to secure the mounting bracket. The short slot acts as a window through which the light beam from the fibre-optic cable is allowed to pass. How far apart the fibre-optic cables (3) need to be mounted depends upon the width of the object to be detected.

Screw the fibre-optic mounting bracket to the chassis at the location shown in the photo. The mounting bracket is secured from underneath with the included screws (3 mm hexagon socket head screws). Secure both mounting brackets.

Note

Make sure that the fibre-optic head is in the correct slot!

Now mount the fibre-optic unit (4) to the mounting plate in the Robotino®. Secure it to the threaded holes provided to this end with the included screws (5).

Inductive analogue sensor

The inductive analogue sensor is simply clamped into the permanently installed retainer. Loosen the knurled screw on the retainer to this end. Push the sensor, plug side up, into the hole provided for this purpose, hold it in the desired position and retighten the knurled screw.

A cable with plug is furnished with the sensor. Plug the cable into the sensor and tighten the knurled nut.

Note

Observe the sensor's working range!

6.13

Connecting the sensors

All of the sensors are connected to the I/O interface via the contact strip, and thus to the controller. Please refer to this documentation and the data sheets for the individual sensors regarding pin allocations for the I/O interface and wiring of the sensors.

6.14 C++ libraries

Different programming options are available to you via C++ to individually modify or expand the functions of the Robotino®.

The CD-ROM provided offers C++ libraries for programming with MS Visual Studio 2005 or higher. These functions enable you to create your own programs to control communication with the Robotino® and the controller of the Robotino® from the PC.

You will find the description of the functions and libraries in the comments within the program.

The PC 104 of the Robotino® contains the same libraries, functions and source data for the Robotino® sample programs in a Linux version. These can be edited with the Linux editor.

Basically two options are available for accessing the functions, libraries and sample programs on the PC 104.

- One option is to access the Linux environment of the PC 104 using a terminal program via a WLAN connection.
- The other option is to connect the keyboard and screen directly to the PC 104 and to access the Linux operating system in this way.

Access with the help of a terminal program

The CD-ROM provided contains the Putty terminal program (putty_0_58.exe), with the following:

1. Establish a WLAN connection to the Robotino®.
2. Start the putty_0_58.exe program. Enter the IP address of the Robotino® in the address line (host name or IP address), and click onto the Open button. The connection with the PC 104 of the Robotino® is established.
3. Enter "Robotino®" as login name at the prompt. Confirm this with the Enter key. You are then requested to re-enter the password. Re-enter "Robotino®" and confirm your entry again.

Access via the system interface of the PC 104

1. Connect your keyboard via one of the two USB ports and a monitor to the VGA interface of the PC 104. The connection is thus established.
2. Enter "robotino" as login name at the prompt. Confirm your entry via the Enter key. You will then be prompted to enter the password. Re-enter "robotino" and confirm your entry again.

In both cases you will then be in the home directory of the user "robotino" of the Linux environment of the PC 104. From here you can access all programs, source codes and libraries.

The directory "examples" contains the sample programs of the Robotino® with their source data and a compiled version.

Example

The program Circle

Change into the directory /robotino/examples/circle and execute the program by entering "./circle" at the prompt. The command starts the already compiled sample program for moving along a circular path (incidentally the same as the one you can start via the display).

By using an editor (e.g. vim), you can now edit the source code and retranslate it using Make.

7. Content of the CD-ROM

7.1 Operating instructions and Data sheets All below mentioned documents are located in the " \Doc\EN" folder on the provided CD-ROM.

Unit	Documents
Motors	MotorGR2042(Diagrams).pdf MotorGR2042(TechData).pdf MotorGR2042(Description).pdf
Accumulators	EN_Powerfit_S3124S(Datasheet).pdf
Distance sensors	EN_Distance_Sensor_gp2d120.pdf
Analogue inductive sensor	EN / DE / ES / FR 678411_Sensor_induktiv_analog_M12.pdf
Incremental encoder	RE30(Data).pdf RE30(Description).pdf
Diffuse sensor	369669_Fibre_optic_device.pdf 369682_Fibre_optic_cable_diffuse.pdf 165327_Fibre_Optic_Device_SOEG_L.pdf 165358_Fibre_Optic_Cable_Diffuse_SOEZ_RT.pdf
Collision sensor	SafetyEdges.pdf
WLAN-Access Point	WAP-0004(Manual).pdf WAP-0004(DataSheet).pdf
Webcam, Windows Help File	Webcam Live Users Guide English.chm
Fuses	Fuses_Robotino®_Datasheet.pdf
Control unit	Kontron_M_MOPSlcdSE_MOPSSE_PSTEM111.pdf

Updates

Current information and amendments to the technical documentation of the Robotino® is available on the Internet at:
<http://www.festo-didactic.com>